USER MANUAL
Bath grab bars 30 cm / 40 cm
Product code: HOLD-TIM-30 / HOLD-TIM-40

1.Characteristic
The bath grab bar is made of durable water-resistant plastic. The surface of the grab bar is corrugated which guarantees a firm grip even in wet conditions. The grab bar can be mounted vertically or horizontally. Mounting hole covers are included.

The set includes: grab bar, washers (2 pcs), wall plug with screw (6 sets).

2. Application
Bath grab bar is ancillary equipment which facilitates tilting to erect position and keeping balance. It is recommended for people with reduced mobility that prevents or hinders the unassisted performance of daily living activities. The grab bar also increases the safety of elderly people using sanitary facilities. It is recommended for use in showers, baths, bathrooms and other humid rooms. For use at home and in hospital.

3. Assembly
Step 1. Before installation, determine the place where the grab bar is to be mounted. While pressing the bracket, mark 6 mounting holes on the wall.
Step 2. Using a impact drill, drill 6 mounting holes at the marked locations.
Step 3. On both sides of the grab bar, between the bar and the wall, place a washer.
Step 4. Using 6 wall plugs with screws, install the grab bar.
NOTE: the type of wall plugs should be adapted to the type of wall (concrete, brick, hollow block, etc.). The included wall plugs are the universal type.

4. Technical data

Product code	Length / Width	Weight

HOLD-TIM-30	30 cm	0,46 kg
HOLD-TIM-40	40 cm	0,56 kg

5. Notes

Timago International Group does not bear responsibility of improper use of the product, failure to comply safety regulations and misuse.

6. Maintenance

- Clean the grab bar regularly using water, a sponge and mild cleaners for plastics such as soap or mild detergents, wipe dry after each use.
- Do not use for cleaning: strong detergents, abrasives, bleach, gasoline or any solvents.
- Check systematically: tightening of screws, condition of grab bar (for cracks)


7. Storage and transport conditions

The product should not be exposed to direct sunlight, water or moisture.

8. Warranty information

All products distributed by our company are covered by a warranty, the terms of which are described in the warranty card available on our website. Please contact the retailer from whom you purchased the product. Please, keep in mind that for warranty purposes, the proof of purchase (receipt or invoice) should be retained.

9. Labels

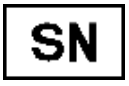	Reference number		Manufacturer
	LOT number		UDI code
	Serial number		Medical device
	Note		Please read the instruction

CE	The manufacturer checked for compliance with essential medical devices requirements.

Our complete line of products is available from our distributors. Learn more at: **www.timago.com**. Thank you for choosing Timago!

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