
USER MANUAL

IV stand with wheels, foldable, 4 hangers

Product code: STAND-TIM

1. Characteristics

The IV stand is made of aluminium. It features four durable hooks and a plastic base with five wheels with brakes. Thanks to its height adjustment mechanism, the stand can be easily tailored to the user's individual needs. It is also foldable, which makes it convenient for transport and storage.

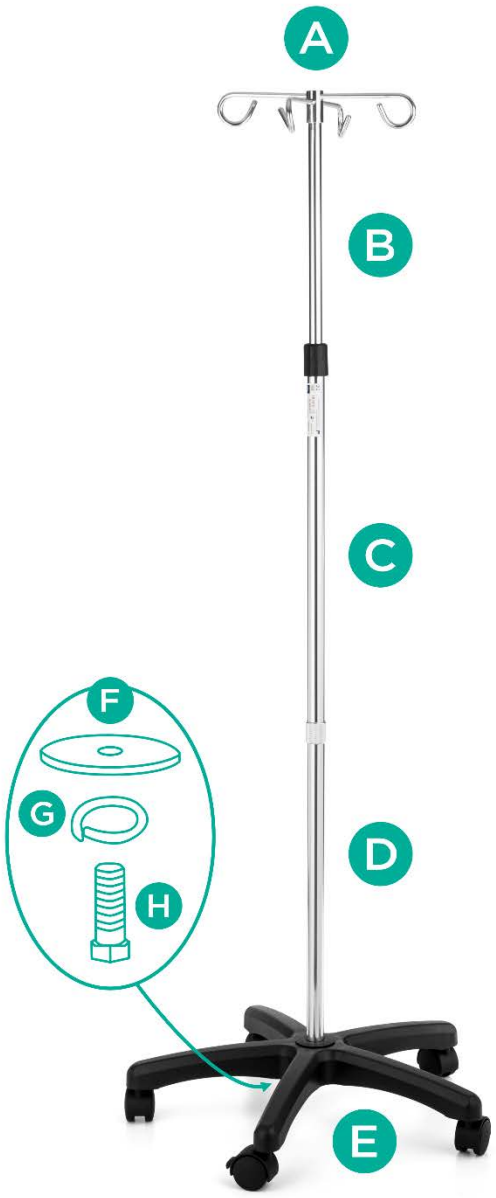Elements of the IV stand:

- **A.** Hooks
- **B.** Top frame tube
- **C.** Middle frame tube
- **D.** Lower frame tube
- **E.** Base with wheels
- **F.** Spacer
- **G.** Spring washer
- **H.** Screw

2. Application

The IV stand with wheels is ancillary equipment designed to facilitate the administration of infusion fluids. It is intended for use both at home and in hospital settings.

3. Assembly and adjustment

Assembly

Step 1. Insert lower tube (D) into the mounting hole located at the base (E).

Step 2. From the bottom of the base (E), fix the spacer (F), the spring washer (G) and tight the screw (H).

Step 3. Inside the bottom tube (D) there's another tube marked with "1" inserted that needs to stay for rest of the assembly process there**.** From the second tube included in set and marked with "1", remove the black securing cap. Then screw together the two tubes marked with the number "1". After connecting two tubes

together, the top tubing (B) will be ready.
Step 4. Overlap the middle tubing (C) on top of the previously prepared top tubing (B) and screw it into the lower tubing (D). The middle tubing (C) is a tube marked with the number "2".
Step 5. Lift the top tubing (B) slightly upward and screw the hooks (A) onto it. Tighten the plastic sleeve on middle tube (C), locking the height of the stand.
The IV stand is ready for use.

HEIGHT ADJUSTMENT
The plastic sleeve in the center of the IV stand is used to adjust the height of the entire IV stand. To adjust, utighten the plastic sleeve, pick the appropriate height and then tighten sleeve afterwards.

BRAKES
The IV stand has 5 wheels with brakes, which allow you to lock the IV stand.

4. Technical data

Height adjustment	119-218 cm
Weight	2,54 kg
Max. capacity of one hanger	2,50 kg
Base diameter	56 cm
Base height	13 cm
Wheels	5x 2'' in diameter

5. Notes
Timago International Group does not bear responsibility of improper use of the product, failure to comply safety regulations and misuse.

6. Maintenance
- Clean the IV stand using water, sponge and mild disinfectants, then dry thoroughly. Do not use for cleaning: strong detergents, abrasives, bleach, gasoline, solvents.
- Check systematically: tightening of mounting elements, condition of hangers with hooks, condition of wheels with brakes.


7. Storage and transport conditions

The product should not be exposed to direct sunlight, water or moisture.

8. Warranty information

All products distributed by our company are covered by a warranty, the terms of which are described in the warranty card available on our website. Please contact the retailer from whom you purchased the product. Please, keep in mind that for warranty purposes, the proof of purchase (receipt or invoice) should be retained.

9. Labels

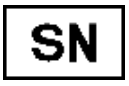	Reference number		Manufacturer
	LOT number		UDI code
	Serial number		Medical device
	Note		Please read the instruction
	The manufacturer checked for compliance with essential medical devices requirements.		

Our complete line of products is available from our distributors. Learn more at: **www.timago.com**. Thank you for choosing Timago!

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