INSTRUCTIONS FOR USE

Patella strap

Product code: TGO-C SR 516

1. Characteristics

Patella strap consists of a strap fastened by Velcro and a silicone pad which affects the patellar ligament with controlled pressure. The insets of the silicone pad perform a micro-massage while movements of knee joint reducing the pain. The strap relieves the patellar ligament and reduces the impact force of quadriceps femoris on tibial attachment.

2. Application

Patella strap is used in medical treatment in order to relief the ligament of the patella.

Indications for use:

- Tears and strains of tibial attachment of quadriceps femoris
- Pathological changes of the patella
- Fractures of the patella
- In non-invasive and preventive treatment of overloads and injuries

Contraindications:

Do not use in case of:

- Skin abrasion
- Eczema
- Wound and other skin injuries
- After application of warming gel or ointment
- Individual contraindications which were determined by the doctor during interview and physical examination

3. Notes

Timago International Group does not bear responsibility of improper

maintenance of the product, failure to comply safety regulations and misuse.

4. Fitting instruction

- Put the strap on with notch directed upwards. The silicone pad should be placed centrally under the knee joint and patella ligament. Wrap the knee joint around with a strap and fasten by Velcro.

Note! Only properly matched size guarantees reduction of pain.

Note! Long-term use may result in reduction of muscle tone and dysfunction of the limb.

5. Technical data

SIZE	UNIVERSAL

6. Maintenance

- Hand wash in warm water. Do not tumble, chlorinate, iron, bleach, dry clean, use tumble dryer.
- Ensure the Velcro is properly fastened in order not to damage the material.
- Dry at room temperature far from intensive heat source.

Note! Maintenance information can be found on the material label.

7. Storage and transport conditions

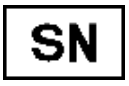
The product should not be exposed to direct sunlight.

8. Warranty information

All products distributed by our company are covered by a warranty, the terms of which are described in the warranty card available on our

website. Please contact the retailer from whom you purchased the product. Please, keep in mind that, for warranty purposes, the proof of purchase (receipt or invoice) should be retained.

9. Labels

	Reference number		Manufacturer
	LOT number		UDI code
	Serial number		Medical device
	Note		Please read the instruction
	The manufacturer checked for compliance with essential medical devices requirements.		

Our complete line of products is available from our distributors. Learn more at: **www.timago.com**. Thank you for choosing Timago!

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