USER MANUAL


Round shower chair, rotating seat, tool free

Product code: TGR-R KP-O 5021

1. Characteristics

The shower chair features an aluminium frame and a rotating seat made of durable PVC. Thanks to its height-adjustable legs, the chair can be adapted to the user's height, while anti-slip tips ensure greater stability and safety. Additionally, the chair is equipped with a Tool-Free (TF) system for easy assembly.

The set includes: seat**,** base tubes (2 pcs), height adjustable legs (4 pcs), mounting screws (5 pcs).

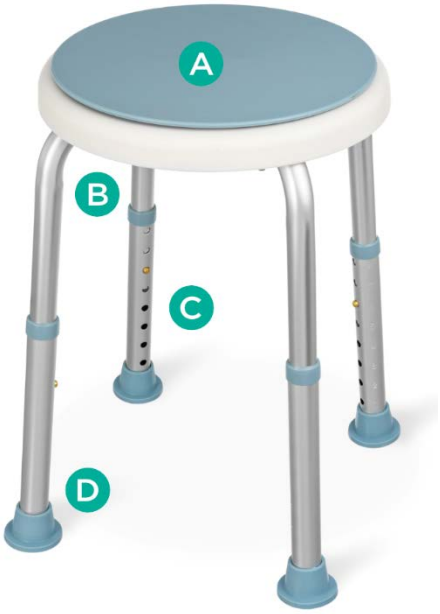Elements of the shower chair:

- **A.** Seat
- **B.** Base tubes
- **C.** Height adjustment – locking pin
- **D.** Anti-sliding tips

2. Application

The shower chair is ancillary equipment which helps the user to perform hygienic activities. The chair is recommended to people with reduced mobility which prevents or hinders the unassisted performance of daily living activities. The chair can be used in the shower, in the bath, in bathrooms and other humid rooms at home and in hospitals.

3. Assembly

Step 1. Place seat upside down on a flat surface.

Step 2. Fix crossed base tubes with seat using four mounting screws with plastic and metal washers. Screw paler wing screw with washers in the middle.

Step 3. Join base tubes with legs.

HEIGHT ADJUSTMENT

- To adjust the height of the chair press the locking pin, select appropriate height of the leg and then lock it again, repeating steps on the other legs. The legs are adjustable in 2.5 cm increments.

Note! All legs must be set on the same height level.

4. Technical data

Height	34,5-52 cm
Seat diameter	32 cm
Weight	1,95 kg
Maximum capacity	110 kg

5. Notes

Timago International Group does not bear responsibility of improper use of the product, failure to comply safety regulations and misuse.

6. Maintenance

- Clean the shower chair regularly using only water, sponge and mild detergents for plastic materials such as soap. Dry thoroughly.
- **Do not use** for cleaning: strong detergents, abrasive agents, bleach, petrol or solvents (it may cause frame damage, visible defects and in consequence - loss of warranty rights).
- At least once every three months check condition of the tubes, locking screws and tightening of screws.
- It is recommended to check rubber tips regularly. In case of worn-out, replace the tips with the new ones.

7. Storage and transport conditions

The product should not be exposed to direct sunlight, water or moisture.

8. Warranty information

All products distributed by our company are covered by a warranty, the terms of which are described in the warranty card available on our website. Please contact the retailer from whom you purchased the product. Please, keep in mind that for warranty purposes, the proof of purchase (receipt or invoice) should be retained.

9. Labels

	Reference number		Manufacturer

LOT	LOT number	UDI	UDI code
SN	Serial number	MD	Medical device
	Note		Please read the instruction
CE	The manufacturer checked for compliance with essential medical devices requirements.		

Our complete line of products is available from our distributors. Learn more at: **www.timago.com**. Thank you for choosing Timago!

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